
ALE-HOP

User Manual / Manual de usuario
Manual d'usuari / Manual do usuário
Mode d'emploi / Benutzerhandbuch
Manuale d'uso / Priručnik s uputama
للتحويل بين الأوضاع

WIRELESS WARM BELT

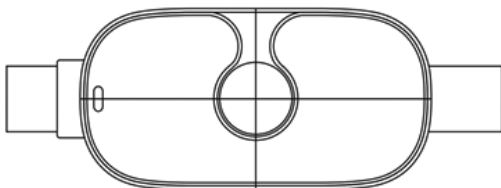Product specifications

Charging power supply: 5V 1A

Electrical power: 5W

Rated voltage: 3.7V


Battery Capacity: 1200mAh


Charging time: 3 hours

Using time: 70-90 mins

Working time: 30 mins auto stop

Product overview

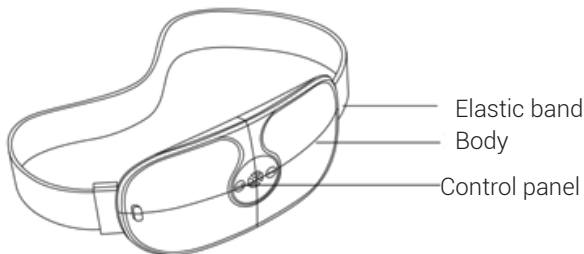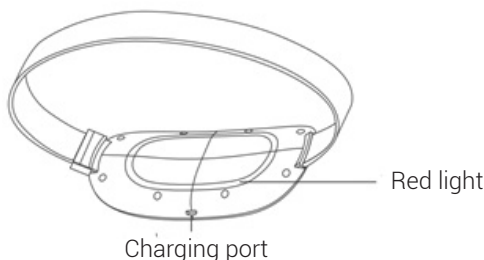Operation and use method:

1. Press and hold the "ON/OFF" button for 2 seconds to switch on the device. The red light will turn on and the heating function will be activated by default.
Press the "HEAT" button once to move to the second heating level.
 2. Press the "HEAT" button again to reach the third level.
Press the "HEAT" button again to reach the third heating level.
 3. Then press the "HEAT" button to turn off the heating. The indicator lights will turn off, and the red light will also turn off.
-
1. Briefly press the "ON/OFF" key to activate the vibration mode.
 2. Press the "ON/OFF" key to enter the second vibration level.
 3. Press the "ON/OFF" button again to enter the third vibration level.
 4. Press the "ON/OFF" button again to deactivate the vibration.
 5. Press and hold the "ON/OFF" button for 2 seconds to switch off the device. The LED indicator will turn off and the device will stop working.

Charging mode

When the heat intensity reduces or the device turns off, it will indicate that the battery is depleted.

1. Connect the charging cable to the charging port and then connect the USB port to a power source (5V).
2. During charging, the charging indicator flashes continuously. When the battery is fully charged, the charging indicator will light up for a long time. It takes about 3 hours to fully charge.
3. If the battery is not charged for a long time, the service life of the rechargeable battery will be shortened.

Caution

1. Carefully read and understand the entire contents of this manual before using this product.
2. Use by pregnant women, children, patients with heart disease, abnormal blood pressure, malignant tumours, cerebrovascular disease, acute patients, or other persons under medical treatment is prohibited.
3. Do not use on sensitive skin.
4. Not for use on patients with skin perception disorders or who are not sensitive to heat.
6. Use the vibration functions during menstruation, as it may cause an increase in menstrual volume. It is recommended to use only the red-light irradiation and heating mode.
7. If you feel uncomfortable or have skin abnormalities during use, please stop use immediately and consult a doctor.
8. Do not use this product in places with high heat, flammability, electromagnetic radiation, and humidity.
9. Do not disassemble the product.
10. It is recommended not to use the product for more than 30 minutes at a time (the product will automatically turn off after 30 minutes). If you need to continue using the product, please turn it on again.
11. Continued use over a long period of time may cause the product to overheat and shorten its life span. If the product overheats, please stop using it and wait for the product to

cool down before using it again.

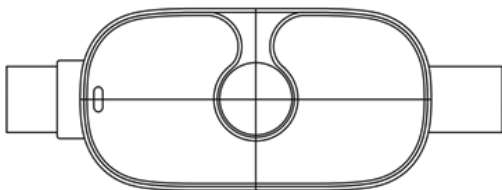
12. When disposing of this product, please do so in accordance with the city's environmental protection regulations.

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Manufactured by | **CLAVE DÉNIA, S.A.**